



Zapper Guide Manual

Read this only if you have done a liver cleanse. Otherwise, this article on parasite cleanse is NOT FOR YOU.

Please read the *Parasite Cleansing* chapter of the book "Cleansing Therapy".

Download it from

<https://www.thetempleofhealing.org/assets/pdf/CleansesFromBook/Parasite.pdf>

and from page 60 to 80 from

<https://www.thetempleofhealing.org/assets/pdf/cleansingTherapy2019.pdf>

Also watch the "*Parasite Cleanse*" YouTube video from by Dr Piyush Saxena

<https://youtu.be/rA-e9EgpQdY>

and for usage by Arjun and Ramesh Pandya Ji at <https://youtu.be/CHxml6HO7xM>.

Clarkia Tincture – 3 times a day

Take 10 drops of Clarkia Tincture (USA made only) three times a day 10 drops in half a cup of water.

<https://www.herbasale.eu/Colon-active-tincture-by-Dr-Hulda-Clark-Parasite-Cleanse-100ml> or sites which sell Clarkia products directly.

Follow suppliers' instructions for usage if they differ from our recommended dosage.

Don't consume made in India parasite tinctures which are sold in online shops like amazon, flipkart etc.

If not available, then take any one of the Homeopathy

- Dr. Reckeweg R56 Worms Drop or
- Wheezal WL42 or
- REPL Dr. Advice No 105 (Worms) or
- Cina Dilution 200 CH

for one week 10 drops three times a day in half a cup of water.

Note: We have used "*Zapper Machine*" to mean "*Laxmi Zapper Pulser*" device.

Consuming Colloidal Silver Water – 3 times a day

Attach Silver Rods with connector and connect this with *Zapper Machine*. Put the Silver rods in a glass full of Water (do not use metal pot/bottle) in a Glass/ plastic



pot/bottle. Press Mode button of *Zapper Machine* 7 times to get Silver Mode. Press Start button once to get it on for 10 min. The colloidal silver water maker has started working now. After 10 minutes the colour of the water will turn milky (for more concentrated and darker Coloidal Silver run it again for 10 min). Medicine is ready.

Switch Off the *Zapper Machine*. Detach silver rods from *Zapper Machine*. Clean them and pack it.

Drink 50 ml of Colloidal Silver water three times a day. For more information about dosage read article on dosage in book "Cleansing Therapy". Colloidal Silver is sensitive to light. Keep it in a coloured glass bottle and store it in a dark place e.g., your refrigerator.

Using Pulsar – 3 times a day

Before Starting:

Make sure that the *Zapper Machine* Power is Off.

Connect the Copper Pipes with *Zapper Machine*.

Rotate the rotary/regulator to the extreme left.

Take two small towels. Soak them in plain salty water.

Instructions for 4 Hz

Switch on the *Zapper Machine*, press Mode button once to get Pulsar **4 Hz**. Start → it.

Wrap the wet towels on the Copper Pipes and hold the Pipes in **both hands**.

Press the Start button in the *Zapper Machine*. Zapping starts.

Slowly rotate the rotary/regulator to the right side, it will start increasing the voltage. Keep rotating to increase current as much as you can bear.

When the Pulsar stops after 10 minutes then switch off the *Zapper Machine*, rotate the rotary/regulator to the extreme left.

Relax for few minutes.



Note: Various variants of using Pulser/ zapper are:

1. Copper Rods wrapped with/without wet towels with/without salt in both hands.
2. One Copper Rod is wrapped with a wet towel in the right hand and another Copper rod is on the sole of the left foot sole. After some time change it to alternate hands and foot. If required, tighten the rods with rope or straps.
3. One Copper Rod is wrapped with a wet towel on the left foot's sole and another Copper rod is on the right foot sole. Thus, covering both legs and pelvic portion.
4. Place the Copper Pipes (covered with a wet towel) below the sole of both feet and tighten them with rope/straps.
5. Look for locations to keep the pipes e.g. armpits and behind the knees.
6. One towel wrapped pipe below waist and the other towel wrapped pipe by hand above the abdomen at various locations. Keep on moving it.

You may take a glass of dilute buttermilk. It helps.

Instructions for 16 Hz

Switch on the *Zapper Machine*, and press Mode button 2 times to get Pulser **16 Hz**. Start → it.

Other instructions remain the same.

Instructions for 100 Hz

Switch on the *Zapper Machine*, and press Mode button 3 times to get Pulser **100 Hz**. Start → it.

Other instructions remain the same.

Using Zapper

Before Starting:

Make sure *Zapper Machine* Power is Off.

Instructions for 2128 Hz:

Switch on the *Zapper Machine*, and press Mode button 4 times to get **2128 Hz**. Start → it.

It will run for 6-7 minutes.



Other instructions remain the same.

Instructions for 10 KHz:

Switch on the *Zapper Machine*, and press Mode button 5 times to get 10 KHz. Start → it.

Other instructions remain the same.

Instructions for 30 KHz:


Switch on the *Zapper Machine*, and press Mode button 6 times to get 30 KHz. Start → it.

It will run for 6-7 minutes.

Other instructions remain the same.

Note: This is one cycle, repeat it three times a day. Do it at least about 3 days.

These frequencies of 4 Hz, 16 Hz and 100 Hz Pulser are square wave 9 v and 2128 Hz, 10 KHz and 30 KHz Zapper are sinusoidal waves. All these 6 frequencies kill different range of parasites. In 3 days, many of them are killed in egg, larva, pupa and adult parasite i.e., all 4 stages of parasites' life.

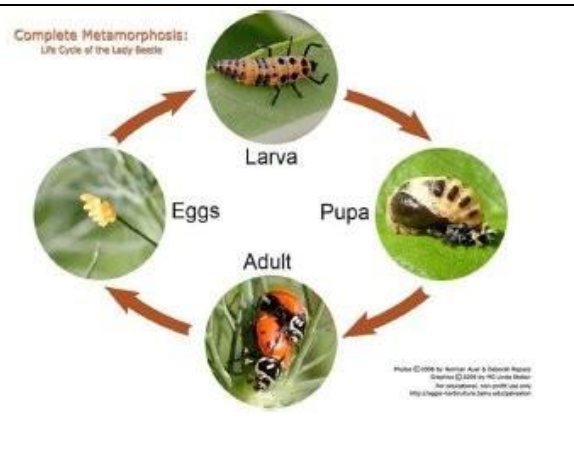

User guide in Hindi is separately available.

This box contains:

1. One Laxmi Zapper Pulser device
2. One Pair of Silver Rods and socket/port. It is good for 5000 glasses of water.
3. One Pair of Copper Pipes with wire and socket/port
4. One Power Cable along with 9V/1A DC Adapter

Buy yourself scrubber, scotch-brite, scrub pad and pitambari for cleaning copper pipes and silver rods.



Few ways to use *Zapper Machine*:



Copper Rods are in Hands



One Rod on Right hand another on Left Sole



One Rod on Left hand another on Right Sole



Both Rods are under the Soles



One Rod is below the back and another is on stomach
(Keep changing the position on stomach)



One Rod is on left side of stomach another is on right side of stomach

<!-- EOF -->